



8-7-3. Movie data (AVI)

You can save movie in specified range in AVI format. Hold down the “Ctrl” key and drag mouse from left to right on waveform to specify time range.

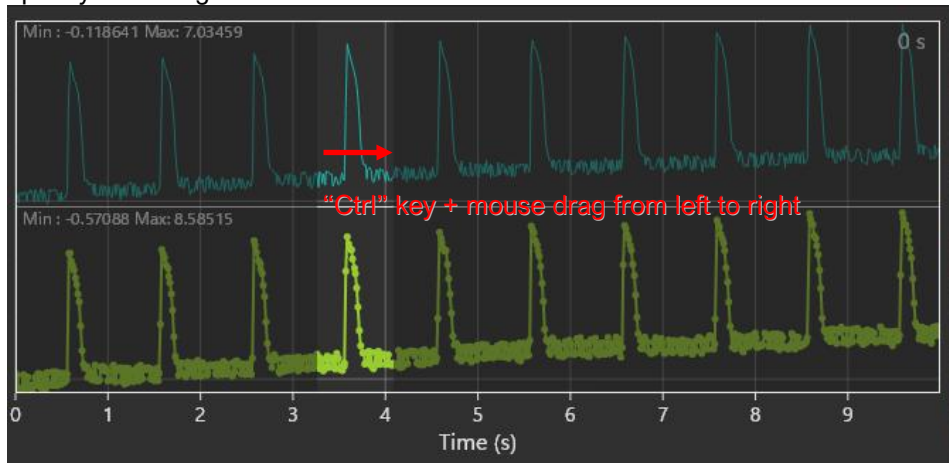

Right-click on each layer image and select [Export video] from the displayed menu.

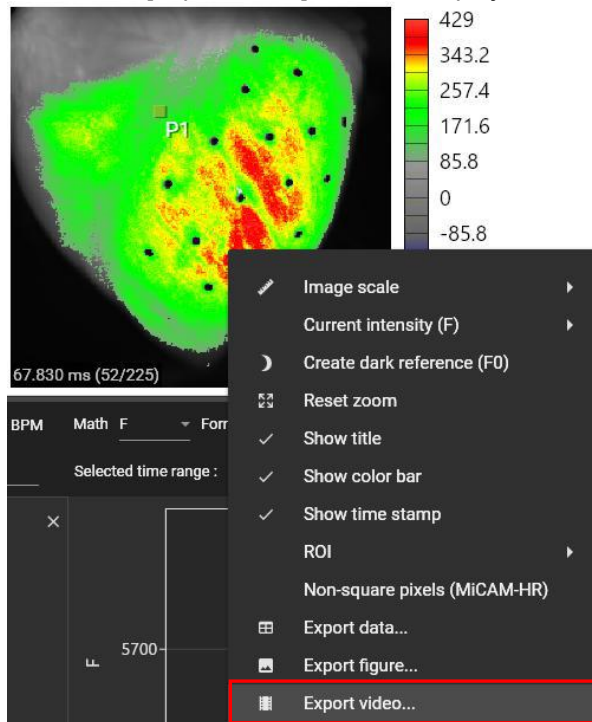

Specify output destination

Display/specify frame range

Specify output movie frame rate

Output movie total time

Check when compressing movie and specify image quality

Output path	C:\Users\T1700forMC05\Desktop\test.avi	
Frame range	16	131
Output frame rate	30	FPS
Output length	00:02:10	
<input type="checkbox"/> Compression	Quality 100	
<input type="checkbox"/> Draw waveforms		

WRITE
CANCEL

Check if waveforms are also displayed

Output movie

Cancel